


Expt. No.: 2

Date: 15.12.19

Experiment	Water Purification Hardness Estimation by EDTA method and its Removal using Ion-exchange Resin
Problem definition	Hardness of water is due to the presence of dissolved calcium and magnesium salts in water. EDTA forms stable complex with hardness causing salts and is used in the removal of scale and sludge forming impurities in industrial boilers.
Methodology	EBT indicator-Metal ion complex is weaker compared to EDTA-metal ion complex. The end point is the color change from wine red (EBT-Metal ion complex) to steel blue (free EBT indicator).
Solution	Estimation of Calcium hardness (in ppm) in the given unknown sample. Understanding the water softening using ion-exchange resins.
Student learning outcomes	Students will learn to a) perform complexometric titration b) understand the efficiency of ion-exchange resins using in water purifiers

Principle:


Ethylenediaminetetraacetic acid (EDTA) forms complexes with a large number of cations including Ca^{2+} and Mg^{2+} depending upon pH of solution. Hence, it is possible to determine the total hardness of water using EDTA solution. EDTA in the form of its sodium salt (H_2Y^{2-}) is commonly used in complexometric titration for estimation of metal ion because pure EDTA (H_4Y) is sparingly soluble in water. EDTA has six binding sites (the four carboxylate groups and the two amino groups) providing six pairs of electrons. The resulting metal-ligand complex, in which EDTA forms a cage-like structure around the metal ion, is very stable at specific pH. All metal-EDTA complexes have a 1:1 stoichiometry. The H_2Y^{2-} form complexes with metal ions as follows.

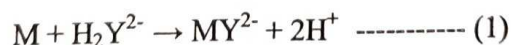

Where, M is Ca^{2+} and Mg^{2+} present in water. Reaction (1) can be carried out quantitatively at pH 10 using Eriochrome Black T (EBT) as indicator. EBT forms a wine-red complex with M^{2+} ions which is relatively less stable than the M^{2+} -EDTA complex. On titration, EDTA first reacts with free M^{2+} ions and then with the metal-EBT indicator complex. The latter gives a colour change from wine-red to steel blue at the equivalence point.

Removal of hardness using ion exchange resins (IER): Ion exchange is a reversible process. When hard water is passed through cation ion-exchange resins packed in a narrow column, Ca^{2+} and Mg^{2+} cations in hard water are exchanged with Na^+ or H^+ ions in the

resins. The exhausted resins are regenerated by passing 10% dil. HCl through the column. A typical example of application is preparation of high-purity water for power engineering, electronic and nuclear industries and in household water purifiers.

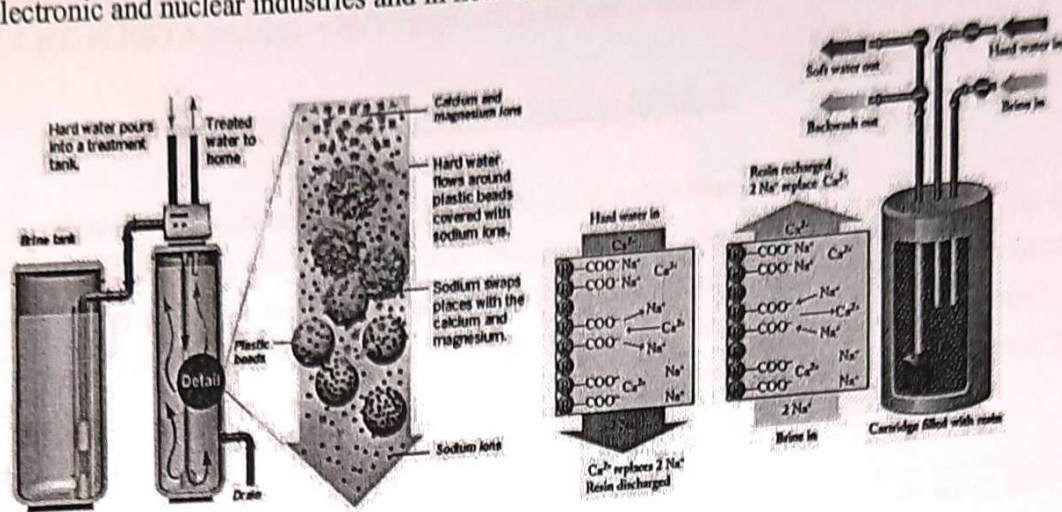

Requirements

Reagents and solutions: Standard hard water (1mg/mL of CaCO₃ equivalents), 0.01 N EDTA solution, EBT indicator, hard water sample, NH₃-NH₄Cl buffer solution and ion exchange resin.

Apparatus: Burette, pipette, conical flask, standard flask burette stand and IER column.

Procedure

Titration-I: Standardization of EDTA

Pipette out 20 mL of the standard hard water containing 1mg/mL of CaCO₃ (1000 ppm) into a clean conical flask. Add one test tube full of ammonia buffer (NH₄OH – NH₄Cl) solution to maintain the pH around 10. Add three drops of Eriochrome Black – T (EBT) indicator and titrate it against the given EDTA solution taken in the burette. The end point is change of colour from wine red to steel blue. Repeat the titration for concordant titer values. Let 'V₁' be the volume of EDTA consumed.

S. No.	Volume of standard hard water (mL)	Burette reading (mL)		Volume of EDTA (V ₁ , mL)
		Initial	Final	
1	20	0	20.2	20.2
2	20	0	20.2	20.2
3				
Concordant titer value				20.2

Calculation:

20 mL of given hard water consumes V_1 mL of EDTA

20 mg of CaCO_3 requires V_1 mL of EDTA for complexation

$\therefore 1 \text{ mL of EDTA requires} = 20/V_1 \text{ mg CaCO}_3 \text{ for complexation}$

This relation will be used in other two titrations

Titration-II: Estimation of total hardness of hard water sample

Pipette out 20 mL of the given sample of hard water into a clean conical flask. Add one test tube full of ammonia buffer ($\text{NH}_4\text{OH} - \text{NH}_4\text{Cl}$) solution and three drops of Eriochrome Black-T (EBT) indicator. Titrate this mixture against standardized EDTA solution taken in the burette. The end point is the change of color from wine red to steel blue. Repeat the titration for concordant titer value. Let ' V_2 ' be the volume of EDTA consumed.

S. No.	Volume of sample hard water (mL)	Burette reading (mL)		Volume of EDTA (V_2 , mL)
		Initial	Final	
1	20	0	16.1	16.1
2				
3				
Concordant titer value				16.1

Calculation:

From Titration 1, we have the following relation:

$\therefore 1 \text{ mL of EDTA requires} = 20/V_1 \text{ mg CaCO}_3 \text{ for complexation}$

From Titration 2,

20 mL of sample hard water consumes = V_2 mL of EDTA.

$= V_2 \times 20/V_1 \text{ mg of CaCO}_3 \text{ eq.}$

$\therefore 1000 \text{ mL of hard water sample consumes} = V_2 \times 20/V_1 \times 1000/20$

$= V_2/V_1 \times 1000 \text{ ppm}$

$\therefore \text{Total hardness of the water sample} = "X" \text{ ppm}$ 797 ppm

$$X = \frac{V_2 \times 1000}{V_1}$$

$$= \frac{16.1 \times 1000}{20.2} = 797$$

Titration-3: Removal of hardness using ion exchange method

Arrange the ion exchange column on to a burette stand and place a clean funnel on top of the column. Pour the hard water sample (around 40 to 50 mL) remaining after the completion of Titration - 2 through the funnel and into the ion exchange column. Place a clean beaker under the column and collect the water passing through the column over a period of 10 minutes. Adjust the valve of the column to match the duration of outflow.

From the water collected through the column, pipette out 20 mL into a clean conical flask and repeat the EDTA titration as carried out above. Note down the volume of EDTA consumed as ' V_3 '.

Calculation:

From Titration 1, we have the following relation:

$\therefore 1 \text{ mL of EDTA requires} = 20/V_1 \text{ mg CaCO}_3 \text{ for complexation}$

From this relation, it can be seen that

20 mL of water sample after softening through the column consumes = V_3 mL of EDTA.

$$= V_3 \times 20/V_1 \text{ mg of CaCO}_3 \text{ eq.}$$

$\therefore 1000 \text{ mL of water sample after softening through the column consumes} =$

$$= V_3 \times 20/V_1 \times 1000/20$$

$$= V_3/V_1 \times 1000 \text{ ppm}$$

$$Y = \frac{V_2}{V_1} \times 1000 = \frac{1.6}{20.2} \times 1000 = 7.9 \text{ ppm}$$

\therefore Residual hardness of the water sample = "Y" ppm

S. No.	Volume of sample hard water (mL)	Burette reading (mL)		Volume of EDTA (V_2 , mL)
		Initial	Final	
1	20	0	1.6	1.6
2				
3				
Concordant titer value				1.6

Result:

- Total hardness of the water sample = "X" ppm = 797 ppm
- Residual hardness in the water sample = "Y" ppm = 7.9 ppm ✓
- Hardness removed through the column = $X - Y$ ppm = 789.1 ppm

Evaluation of Result:

Sample number	Experimental value	Actual Value	Percentage of error	Marks awarded
				10/10.

19BCE2074

ok
19/12/19.